

Continuous freeze-drying and its relevance to the pharma/biotech industry

Original

Continuous freeze-drying and its relevance to the pharma/biotech industry / Pisano, Roberto; Capozzi, Luigi C.; Trout, Bernhard. - ELETTRONICO. - (2017), pp. 1-26. (Intervento presentato al convegno Integrated Continuous Biomanufacturing III tenutosi a Cascais (Portugal) nel September 17-21).

Availability:

This version is available at: 11583/2701006 since: 2018-02-22T09:51:17Z

Publisher:

ECI Symposium Series

Published

DOI:

Terms of use:

openAccess

This article is made available under terms and conditions as specified in the corresponding bibliographic description in the repository

Publisher copyright

(Article begins on next page)

9-19-2017

Continuous freeze-drying and its relevance to the pharma/biotech industry

Roberto Pisano

Politecnico di Torino, roberto.pisano@polito.it

Luigi C. Capozzi

Politecnico di Torino

Bernhardt Trout

Massachusetts Institute of Technology

Follow this and additional works at: http://dc.engconfintl.org/biomanufact_iii




Part of the [Engineering Commons](#)

Recommended Citation

Roberto Pisano, Luigi C. Capozzi, and Bernhardt Trout, "Continuous freeze-drying and its relevance to the pharma/biotech industry" in "Integrated Continuous Biomanufacturing III", Suzanne Farid, University College London, United Kingdom Chetan Goudar, Amgen, USA Paula Alves, IBET, Portugal Veena Warikoo, Axcella Health, Inc., USA Eds, ECI Symposium Series, (2017).
http://dc.engconfintl.org/biomanufact_iii/70

This Abstract and Presentation is brought to you for free and open access by the Proceedings at ECI Digital Archives. It has been accepted for inclusion in Integrated Continuous Biomanufacturing III by an authorized administrator of ECI Digital Archives. For more information, please contact franco@bepress.com.

Continuous freeze-drying and its relevance to the pharma/biotech industry

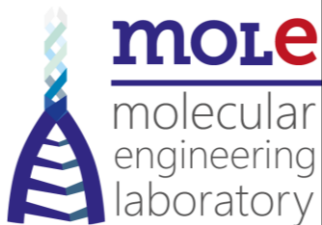

Roberto Pisano

Department of Applied Science and Technology
POLITECNICO di TORINO

Introduction

Freeze-drying /Lyophilization is

a process where water, or another solvent, is removed from a frozen solution at low temperature and pressure via sublimation

A freeze-drying cycle encompasses three steps

- Freezing
- Primary drying
- Secondary drying

Almost 50% of biopharmaceuticals listed by FDA and EMA is lyophilized, proving that freeze-drying is the preferred way to stabilize large molecules that are not stable in liquid, despite its high energy consumptions and long processing time.

Introduction

Freeze-drying of pharmaceuticals is performed batch-wise

- Long and expensive process

vial-to-vial
heterogeneity

- Heat and mass transfer is not uniform within the batch of vials

- Heterogeneity in freezing behavior


- Heterogeneity in drying behavior

- Poor control of product quality



Examples of lyophilized samples
belonging to the same lot of production


Background & Problem statement

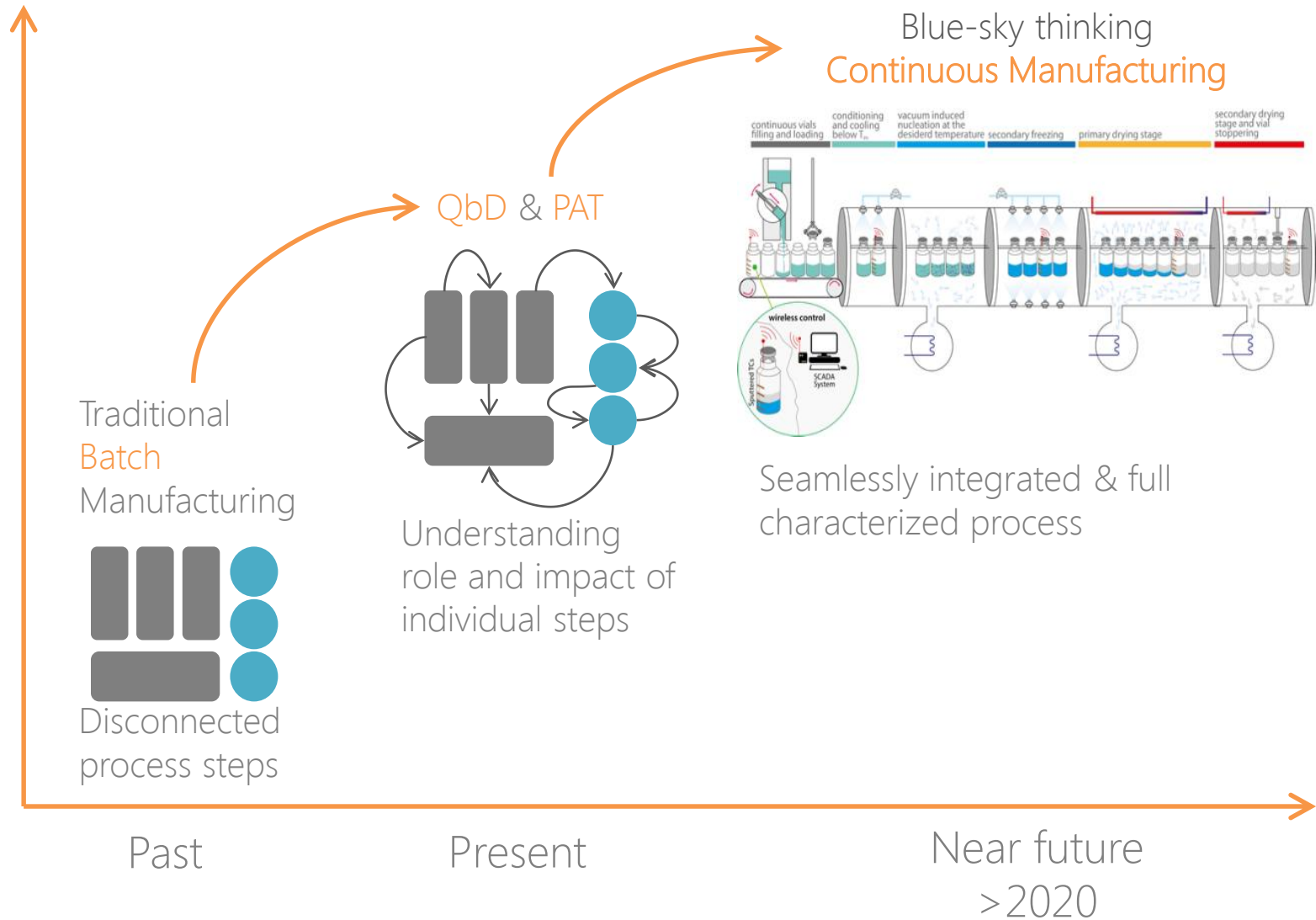

Drawbacks of batch freeze-drying

Heterogeneity in freezing behavior ...

- ▶ temperature of nucleation is not uniform within the batch of vials, but is stochastically distributed,

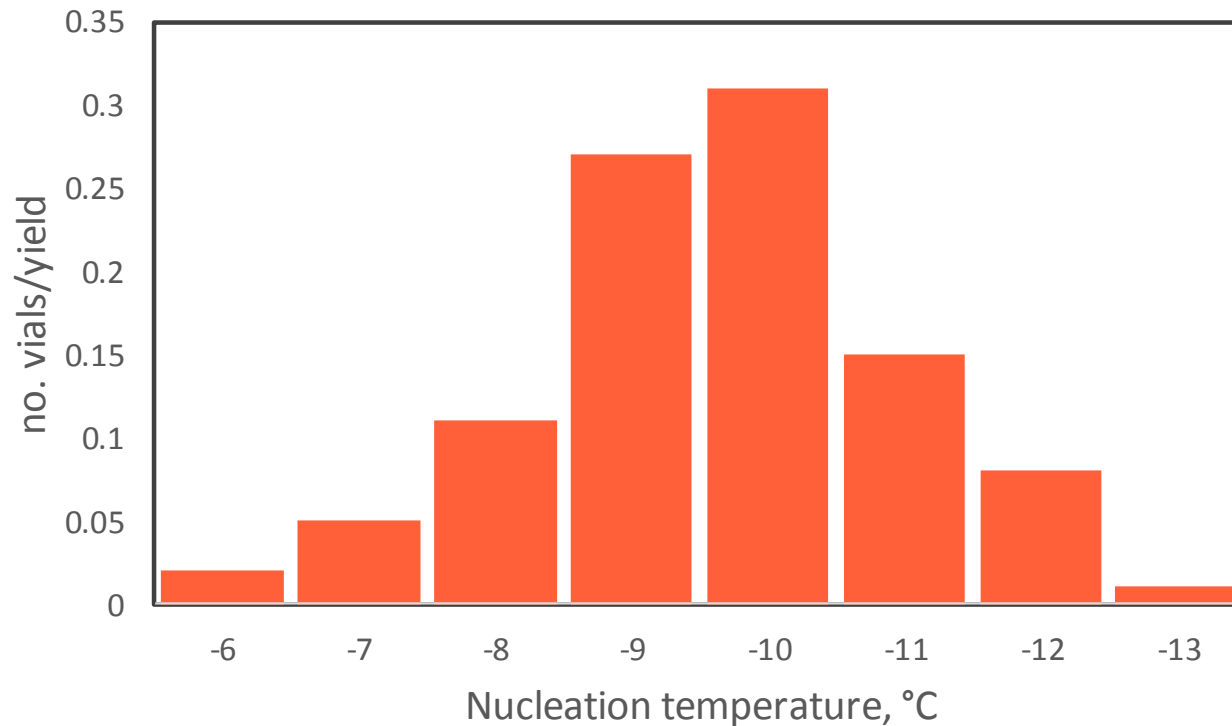


Distribution of the nucleation temperature as observed
in a batch freeze-drying cycle

Drawbacks of batch freeze-drying


Heterogeneity in freezing behavior ...

- ▶ temperature of nucleation is not uniform within the batch of vials, but is stochastically distributed
- ▶ ice structure and, hence, cake morphology changes from vial to vial

$T_n = -10\text{ °C}$

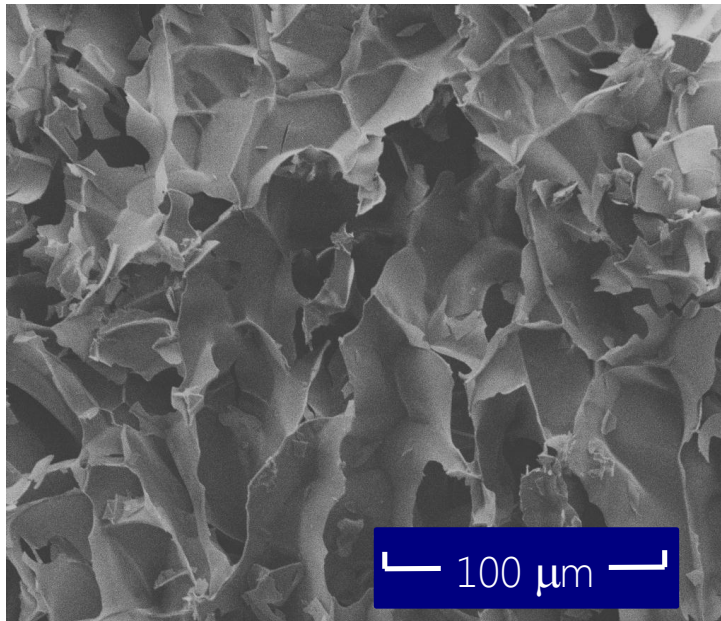

$T_n = -15\text{ °C}$

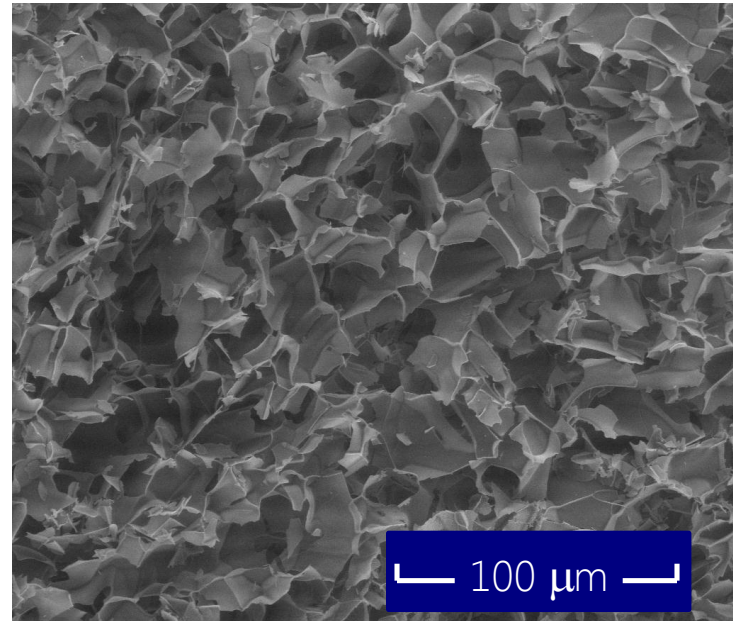

SEM micrographs of mannitol 5% as produced by batch freeze-drying

Drawbacks of batch freeze-drying

Heterogeneity in freezing behavior ...



- ▶ temperature of nucleation is not uniform within the batch of vials, but is stochastically distributed
- ▶ ice structure and, hence, cake morphology changes from vial to vial
- ▶ both primary and secondary drying behavior change from vial to vial
- ▶ vial-to-vial variations in polymorphs composition
- ▶ large distributions in residual moisture and potentially in API activity/stability

Continuous freeze-drying might be beneficial to ...

- ▶ achieve a narrow distribution in nucleation temperature
- ▶ make the frozen product morphology more uniform
- ▶ make drying behavior more uniform among the vials of the batch
- ▶ reduce vial-to-vial heterogeneity

Drawbacks of batch freeze-drying


Heterogeneity in heat transfer ...

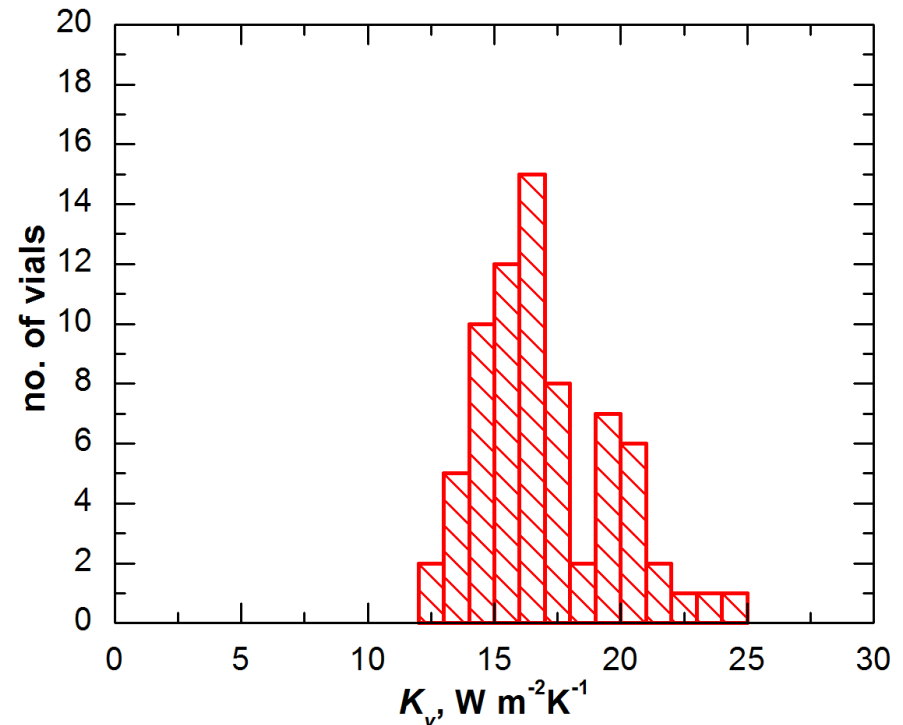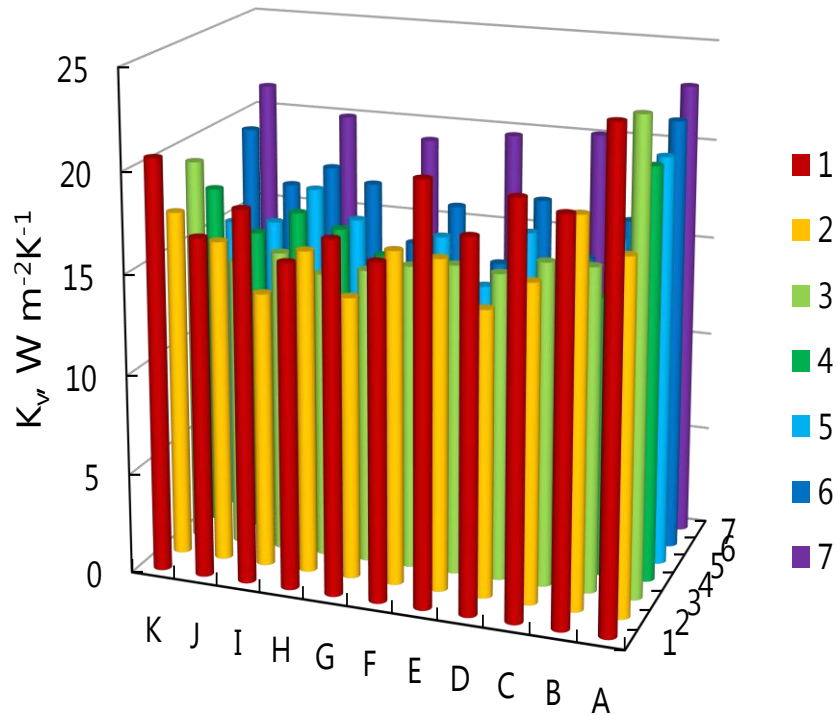


Spatial and statistical distribution of the heat transfer coefficient, between shelf and container, within a batch of vials. Data refer to primary drying, 10 Pa as chamber pressure

Drawbacks of batch freeze-drying

Heterogeneity in heat transfer ...

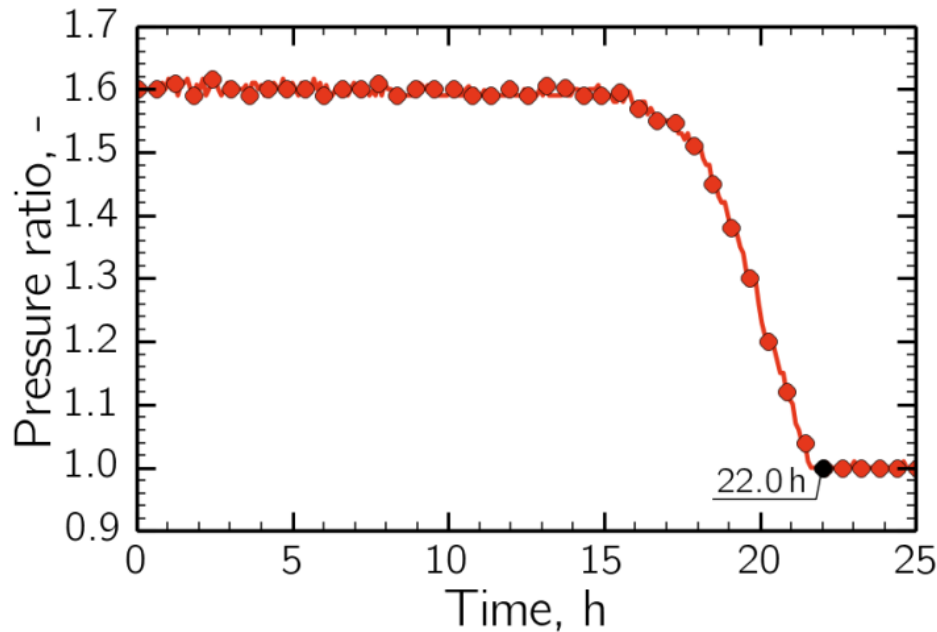


Evolution of pressure ratio as observed
in a batch freeze-dryer

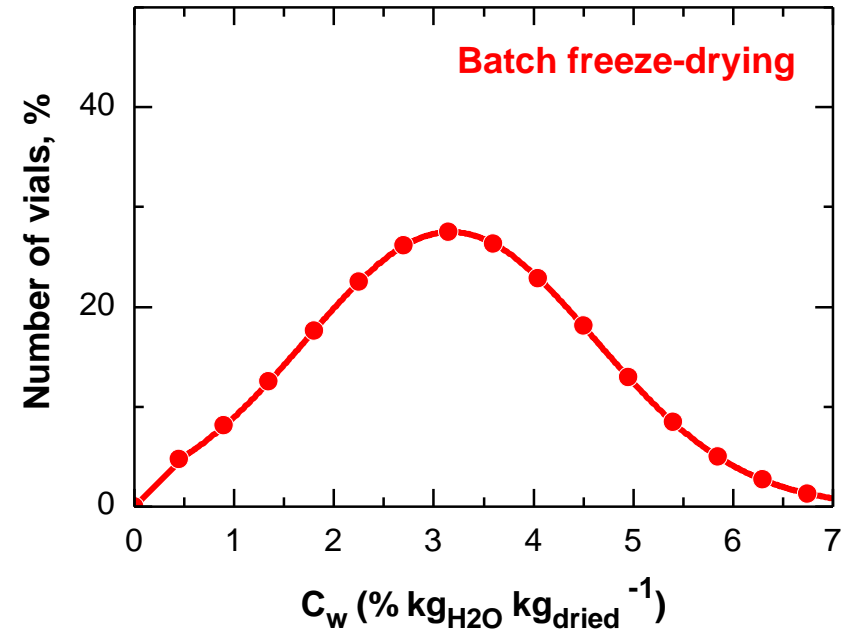

Statistical distribution of the residual
moisture within the lyophilized
samples at the end of primary drying

Drawbacks of batch freeze-drying

- ▶ Variations in product morphology due to freezing
- ▶ Variations in the residual moisture at the end of primary drying
- ▶ Variations in the residual moisture at the end of secondary drying

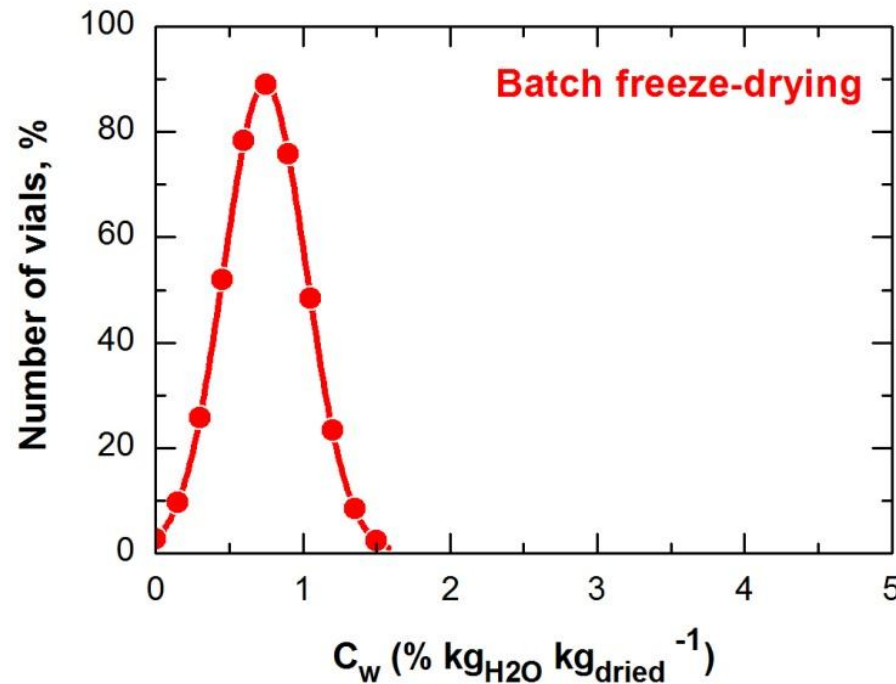


Statistical distribution of residual moisture within the lyophilized samples as observed at the end of secondary drying

The extent of heterogeneity in freezing and drying behavior is equipment-specific.

A cycle developed in a laboratory freeze-dryer cannot be transferred without modifications to the production unit → **scale up**

A new concept for the continuous freeze-drying of unit doses

OBJECTIVE: development of a continuous freeze-dryer that produces a final product having similar properties and structures to that obtained by a conventional batch unit.

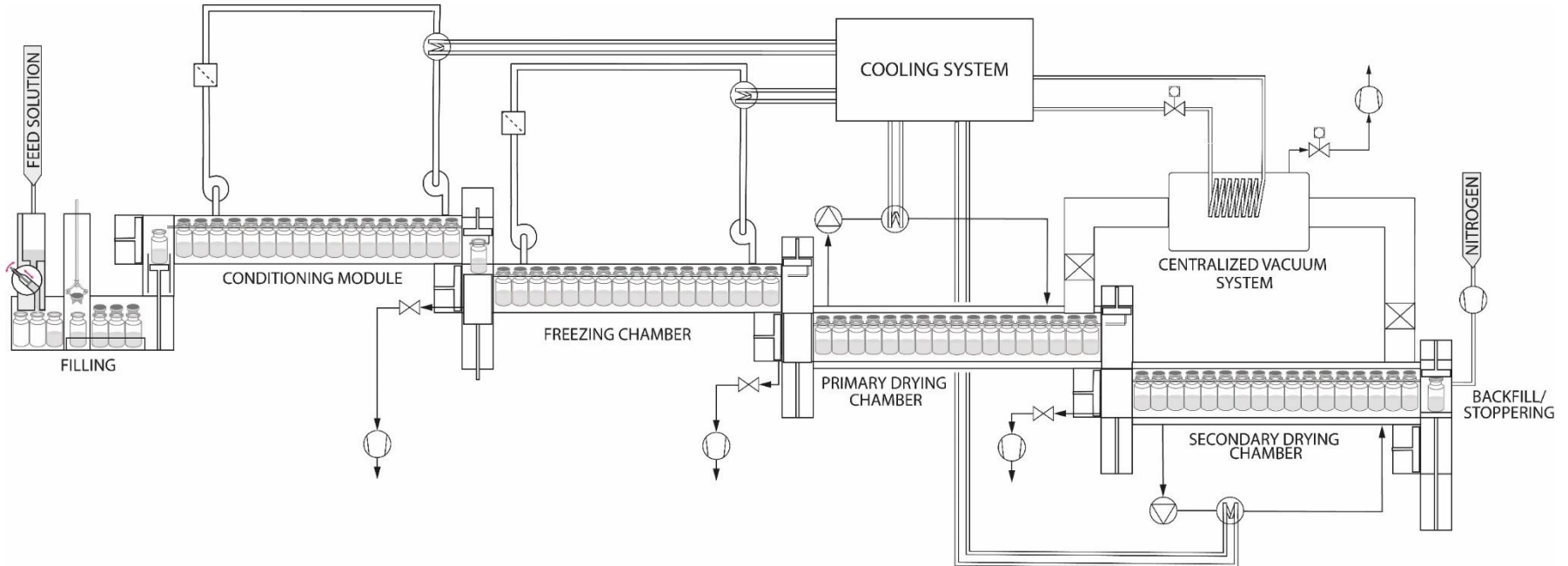

FILLING

CONDITION-
ING

VACUUM
CHAMBER

PRIMARY
DRYING
MODULE


SECONDARY
DRYING
MODULE

BACK/STOP
PERING


A new concept for the continuous freeze-drying of unit doses

Filling and Loading

Conditioning module
Nucleation module
Freezing module
Primary drying module
Secondary drying module

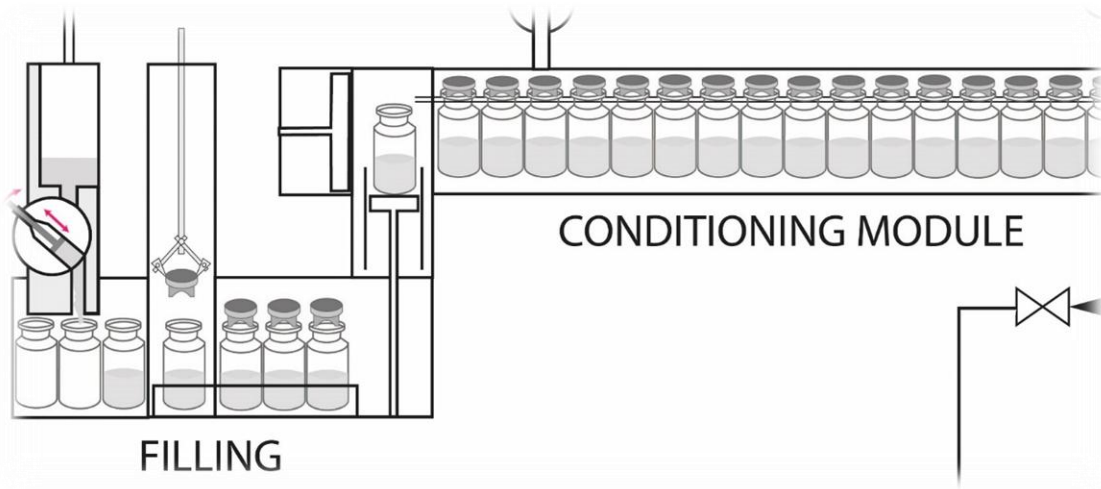

Moving of vials

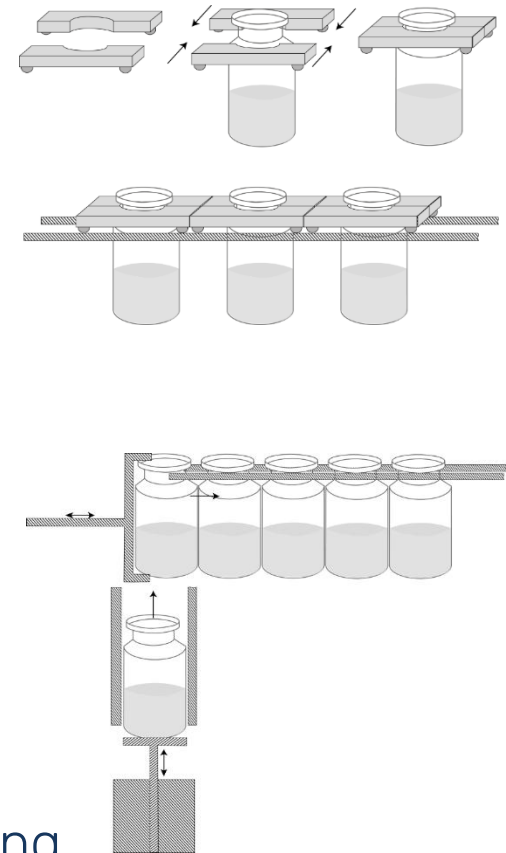

The continuous flow of vials is achieved by suspending the vials over a track → **uniformity in heat transfer**

A new concept for the continuous freeze-drying of unit doses

Filling and Loading


Conditioning module

Nucleation module

Freezing module

Primary drying module

Secondary drying module

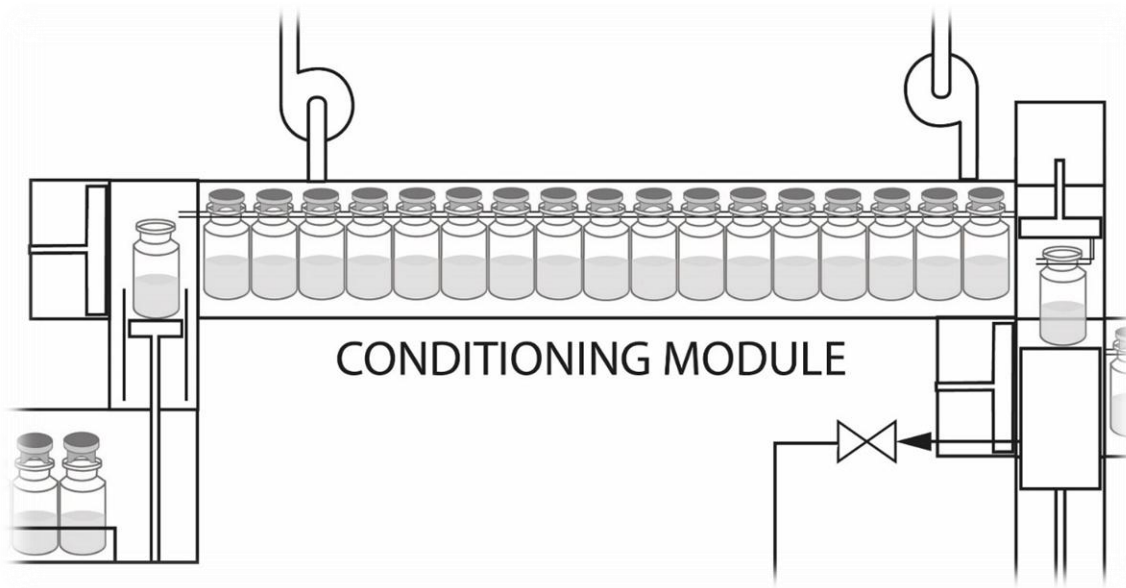


A new concept for the continuous freeze-drying of unit doses



Filling and Loading
Conditioning module

Nucleation module

Freezing module
Primary drying module
Secondary drying module

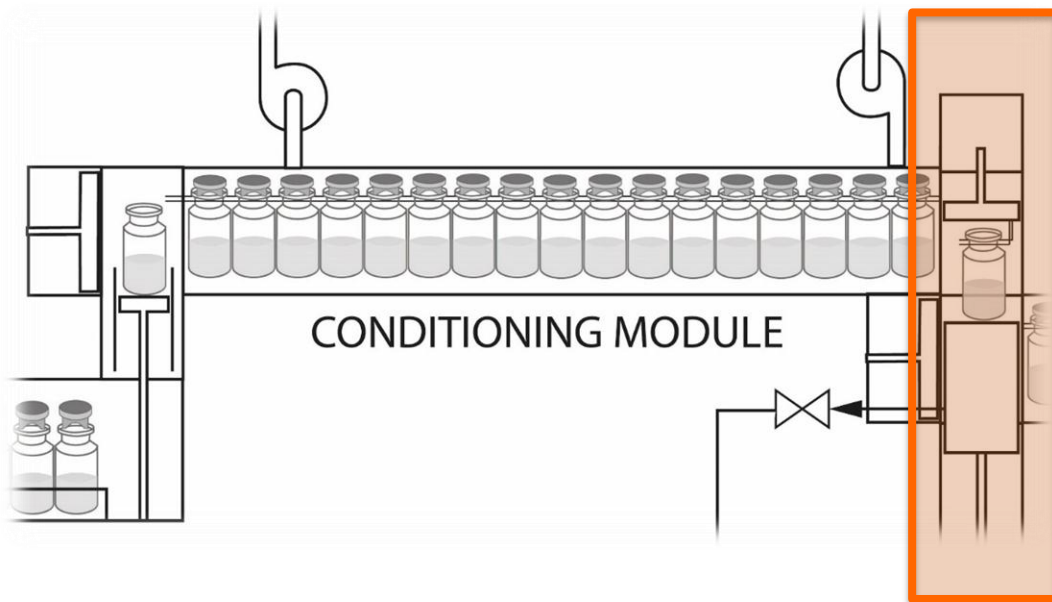

Example of nucleation chambers

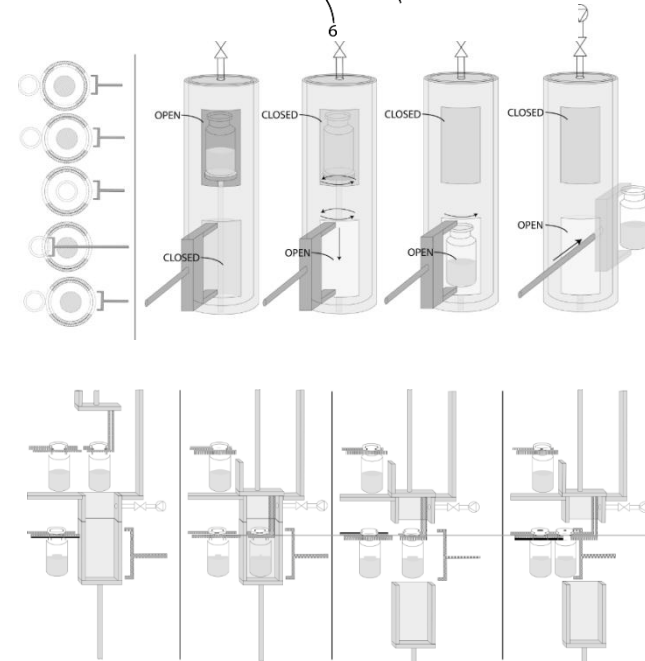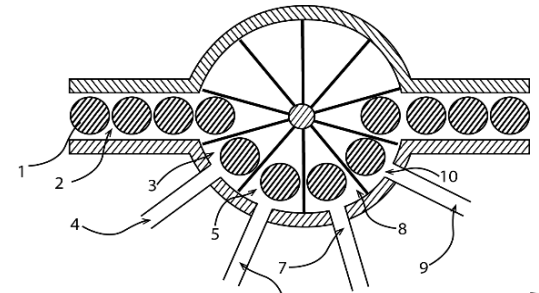


A new concept for the continuous freeze-drying of unit doses

Filling and Loading
Conditioning module
Nucleation module

Freezing module

Primary drying module
Secondary drying module

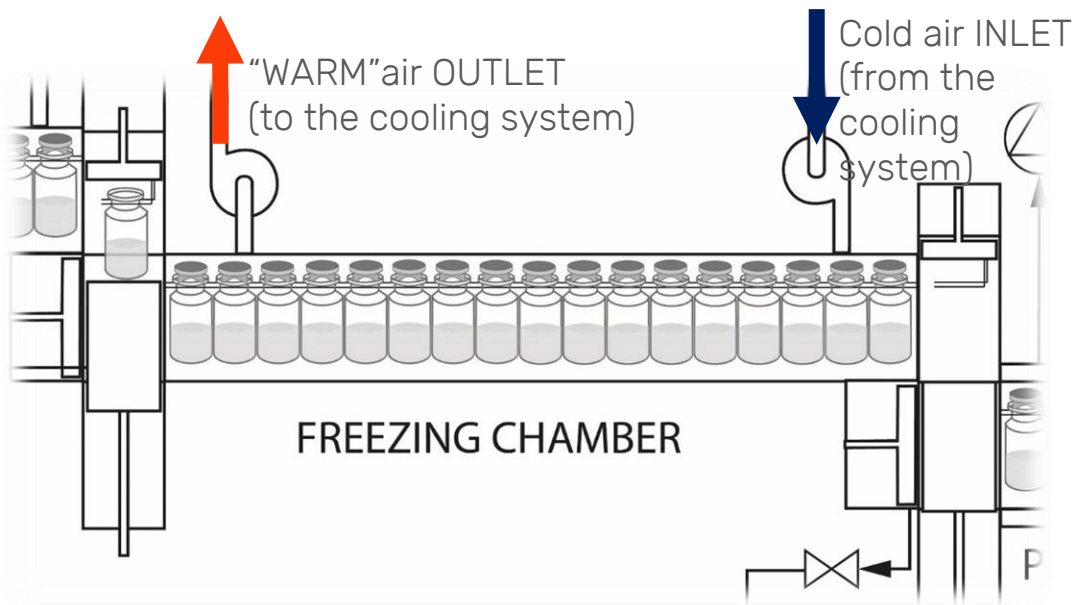


- The nucleated solution is further cooled by forced convection until its complete solidification.
- The external surface of the vessel is equally flushed by the cryogenic gas.
- Different freezing protocols can be performed modulating temperature and velocity of cryogenic gas.

A new concept for the continuous freeze-drying of unit doses

Filling and Loading
Conditioning module
Nucleation module
Freezing module

Primary drying module

Secondary drying module

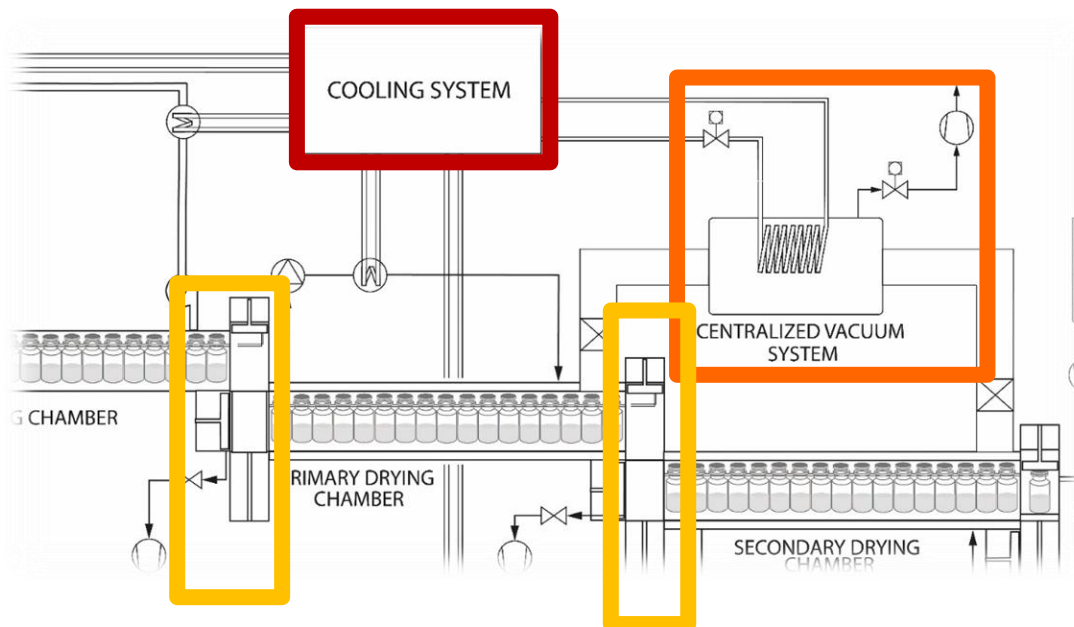

In the primary drying module ...

- Vials are exposed to low temperature and pressure
- Heat is transferred by radiation from temperature-controlled surfaces

- Vacuum system (condenser + vacuum pump)
- Cooling/heating system
- Sluice-gate/load-lock


A new concept for the continuous freeze-drying of unit doses

Filling and Loading
Conditioning module
Nucleation module
Freezing module
Primary drying

Secondary drying module

In the secondary drying module ...

- Vials are exposed to high temperature and low pressure so as to promote desorption of bounded water

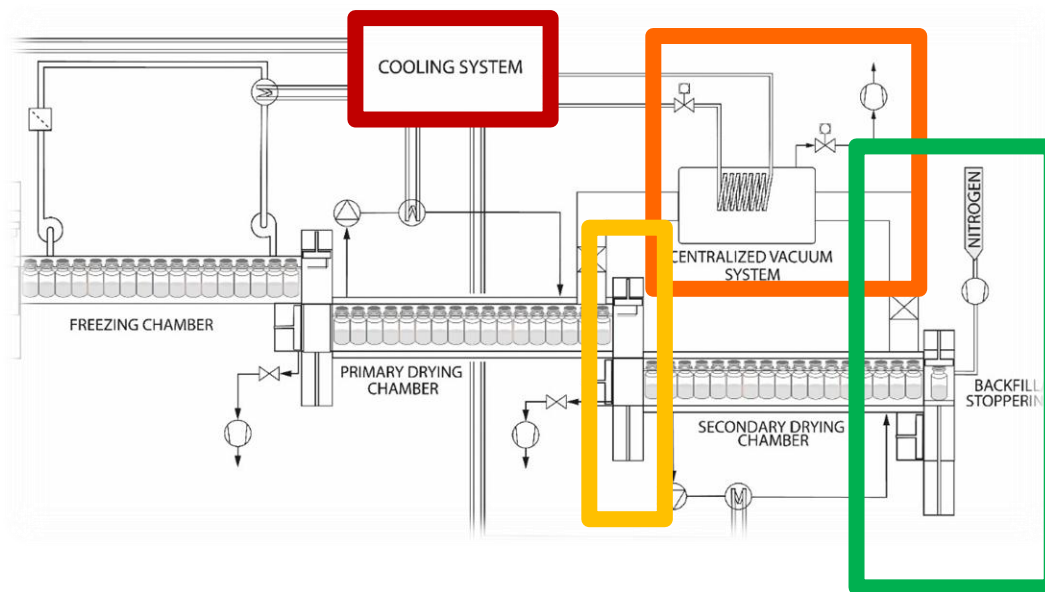


- Vacuum system
(condenser + vacuum pump)
- Cooling/heating system
- Sluice-gate/load-lock
- Stoppering/sealing

A new concept for the continuous freeze-drying of unit doses

Flexibility & Modularity


- The various modules can be combined to make the system more flexible and treating products from different upstream feeds.

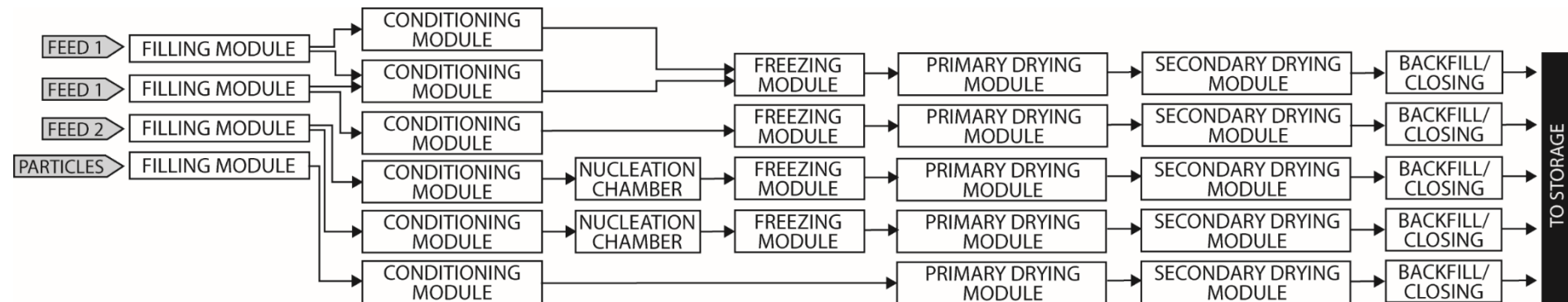

Experimental results

Product morphology

- More precise control of freezing conditions
- Larger pores and hence smaller resistance to mass transfer during primary drying



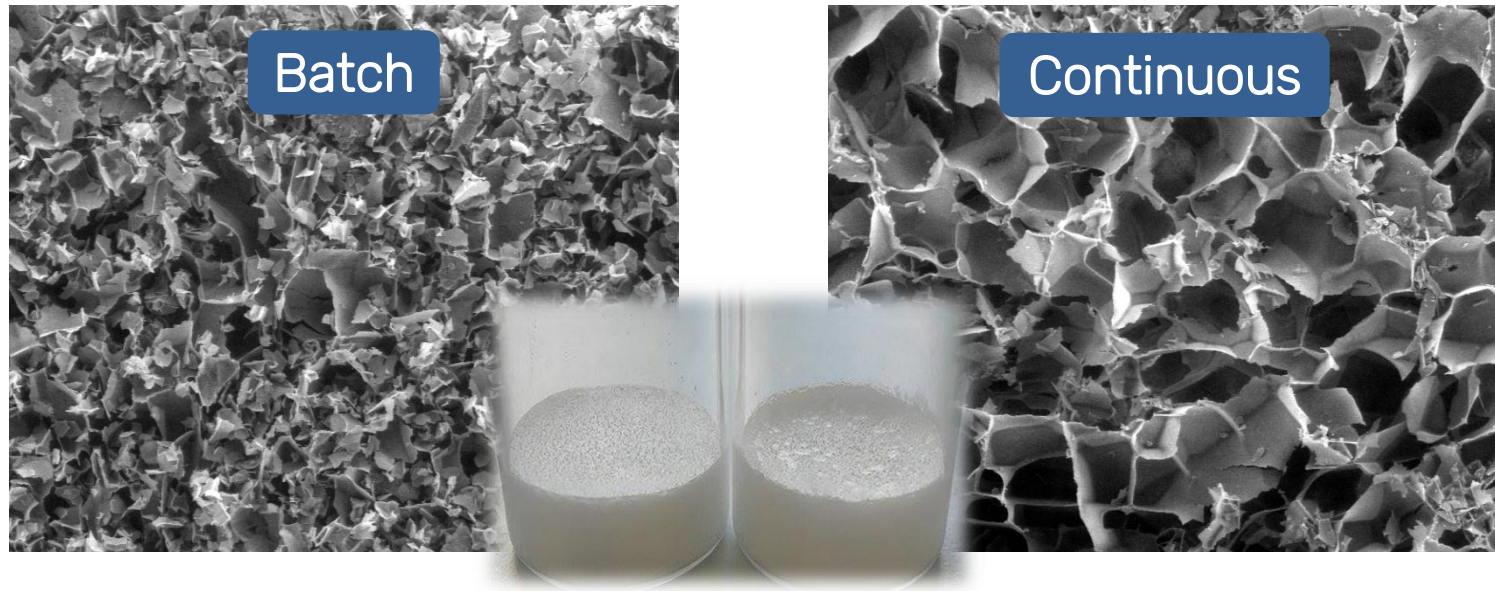


SEM images of lyophilized mannitol samples produced on constant drying conditions. Images refer to the same enlargement



Experimental results

Product morphology

- More precise control of freezing conditions
- Larger pores and hence smaller resistance to mass transfer during primary drying
- Intra-vial heterogeneity is less evident

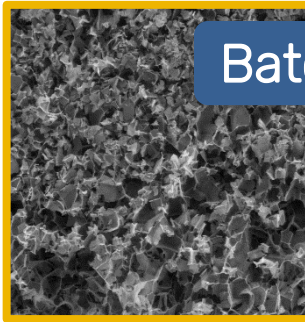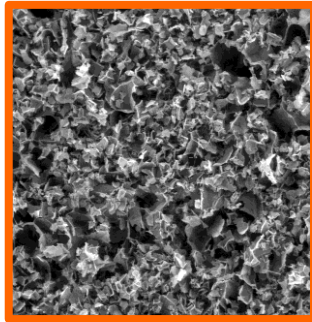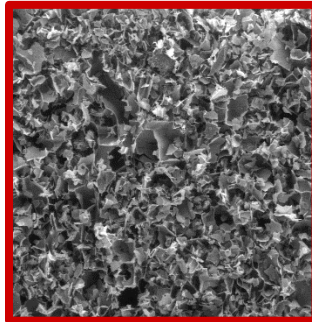

Batch

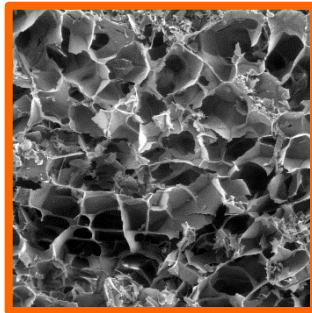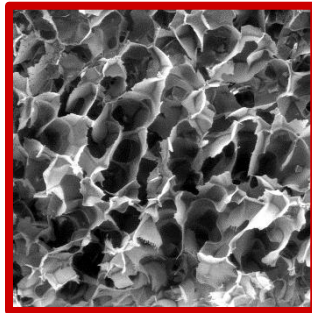

Continuous



SEM images of lyophilized mannitol samples

Experimental results

Process performances


- Larger pores and hence shorter primary drying
- Breaks of a typical batch production can be 20% to 50% of the total cycle time
- The overall cycle time is up to **5 times shorter**

	Loading	Leak test	Freezing	Primary drying	Soak time	Secondary drying	Closing	Unloading	Defrost/CIP /SIP/H2O2
Batch	✓ 5 h	✓ 2-3 h	✓ 6 h	✓ LONG	✓ 6 h	✓ SHORT	✓ 1 h	✓ 6 h	✓ 6 h
Continuous			✓ < 1 h	✓ SHORTER		✓ SHORT			

Experimental results

Process performances

- Larger pores and hence shorter primary drying
- Breaks of a typical batch production can be 20% to 50% of the total cycle time
- The overall cycle time is up to 5 times shorter
- Distribution of the residual moisture at the end of drying is more uniform

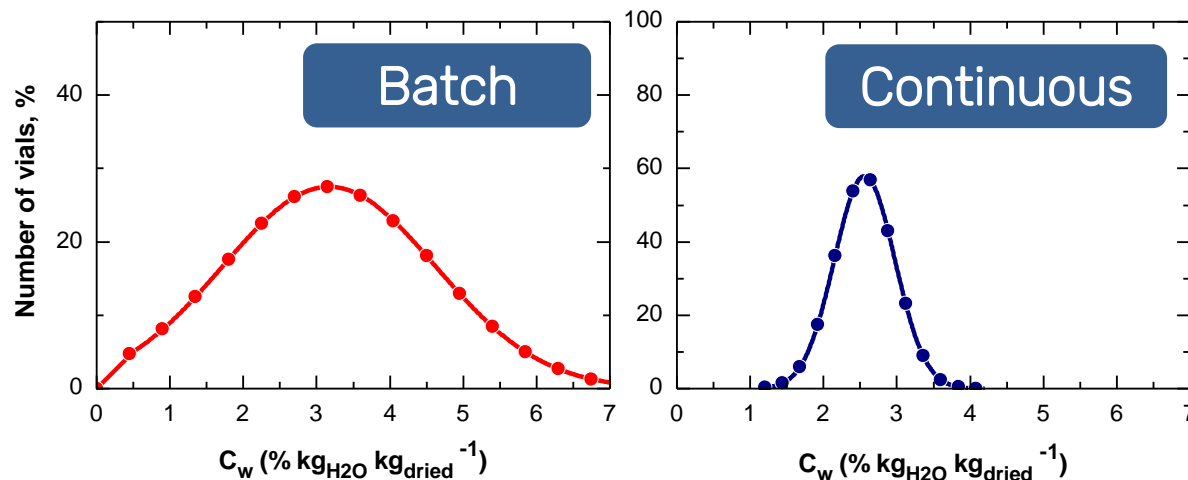


Distribution of the final residual moisture for a sucrose-based formulation (250 vials)

Experimental results

Case study #1 – 100,000 vials/week

The equipment volume is approx. 15 times smaller

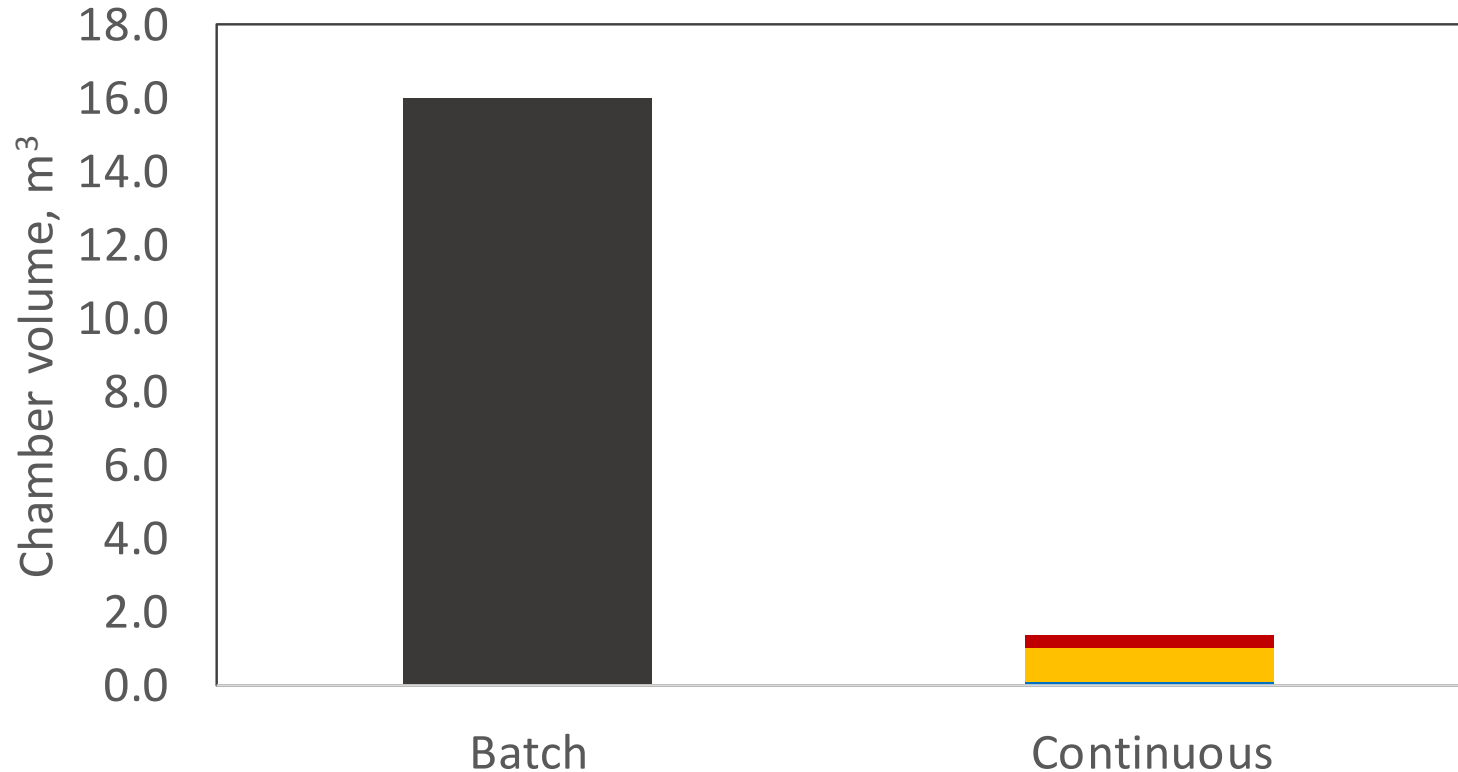

Conclusions

Reduce the risk of product contamination

- No manual handling, increased safety
- The processing time is shorter

Modular and smaller equipment and facilities

- More flexible operation
- Reduced inventory
- Lower capital costs, less work-in-progress materials

Eliminate scale-up from lab to production units

Process flexibility

- Bulk vs. particle-based material
- Yield can be adjusted on market request

Improve product quality

- Uniformity of the lot of production
- In-line control of product quality

Acknowledgements



Politecnico di Torino

Department of Applied Science and Technology

- Andrea Arsiccio
- Fiora Artusio
- Marco Bazzano
- Luigi C. Capozzi



Massachusetts Institute of Technology

Department of Chemical Engineering

Prof. Bernhardt Trout

Thanks for
your attention!

