

S. aureus Colonized Human Skin Equivalent - In Vitro Bio-evaluation Tool for Antibacterial Polymers

Original

S. aureus Colonized Human Skin Equivalent - In Vitro Bio-evaluation Tool for Antibacterial Polymers / Idrees, Ayesha; Viebahn, Richard; Ciardelli, Gianluca; Chiono, Valeria; Salber, Jochen. - (2018). (Intervento presentato al convegno Drug-free antibacterial technology for medical applications tenutosi a The Hauser Forum, Cambridge nel 14 December, 2018).

Availability:

This version is available at: 11583/2722662 since: 2019-03-06T15:06:25Z

Publisher:

Not applicable

Published

DOI:

Terms of use:

openAccess

This article is made available under terms and conditions as specified in the corresponding bibliographic description in the repository

Publisher copyright

(Article begins on next page)

S. aureus Colonized Human Skin Equivalent - *In Vitro* Bio-evaluation Tool for Antibacterial Polymers: ESR14

Ayesha Idrees^{1, 2}, Richard Viebahn², Gianluca Ciardelli¹, Valeria Chiono¹, Jochen Salber²


¹Department of Mechanical and Aerospace Engineering (DIMEAS), Politecnico di Torino, Italy

² UMC Knappschaftskrankenhaus GmbH, Clinic of Surgery, Hospital of the Ruhr-University, Bochum, Germany

A bacterial colonized human skin equivalent (c-HSE)





The aim of this study was the development of a human skin model and human skin wound infection model for the bio-evaluation of antimicrobial biomaterials intended for wound healing purposes. The three-dimensional *in vitro* models will represent advanced and complex systems to perform more reliable preclinical studies. These models will be employed for *in vitro* screening of both antibacterial activity as well as cytocompatibility of new biomaterials and wound dressings to optimize their *in vivo* performances. In the study, the 3D systems were developed and their structure was characterized. The antibacterial activity and cytotoxicity of a model wound dressing releasing Ag⁺ was analyzed in the models.

3D Dermal Fibroblast Model

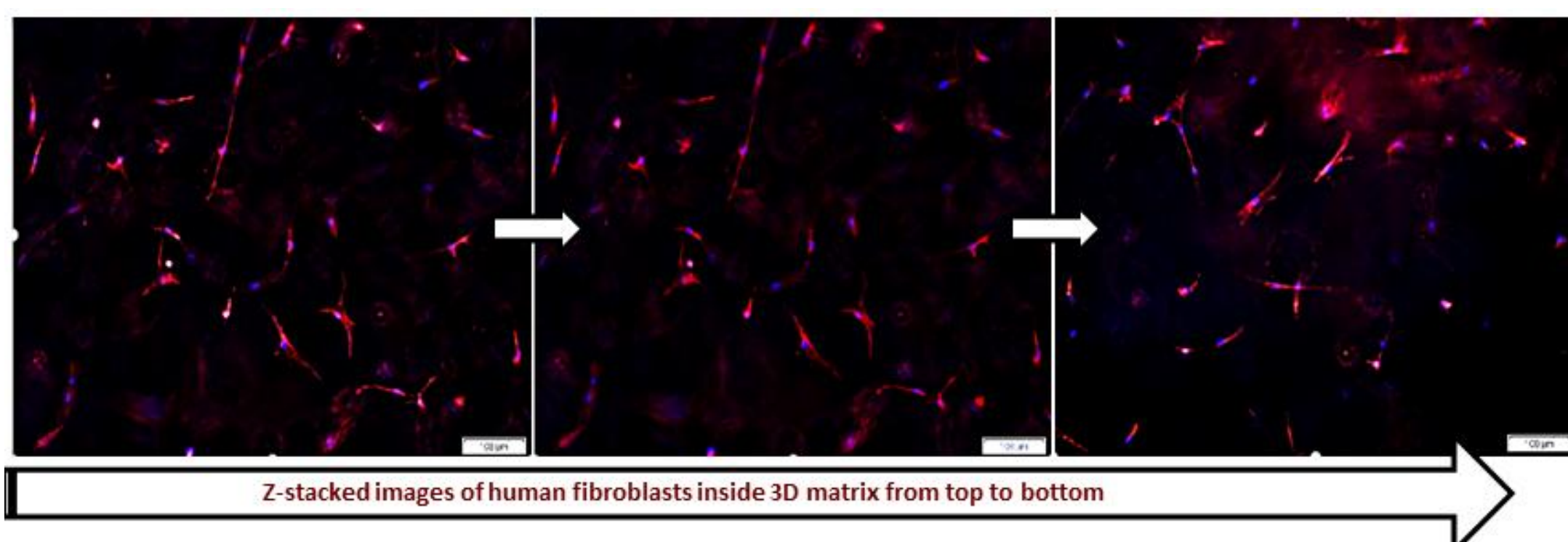

Optimizing the dermal part of human skin: Z-stacked imaging revealed the filopodia like morphology and a uniform distribution of human fibroblasts at different planes inside a Col-I matrix. Fluorescent microscopic images show cell nuclei stained with DAPI and cytoskeletal F-actin stained with Phalloidin. Scale bar=100 µm

Immunohistochemistry of the 3D Human Skin Model

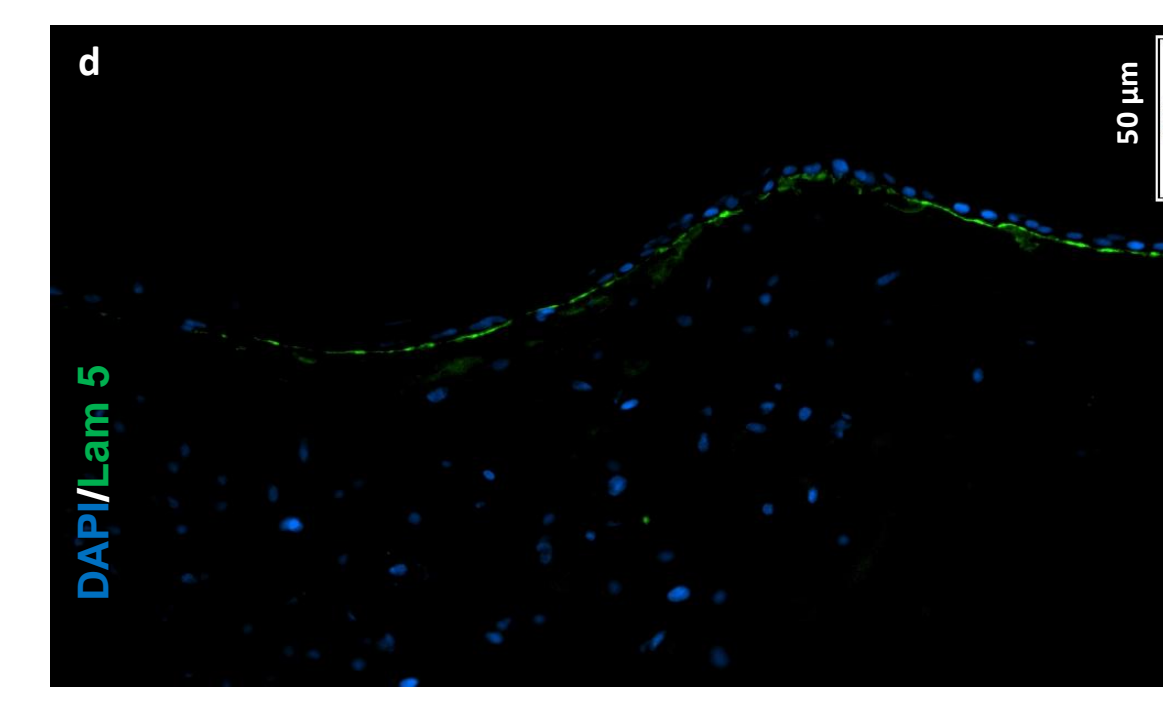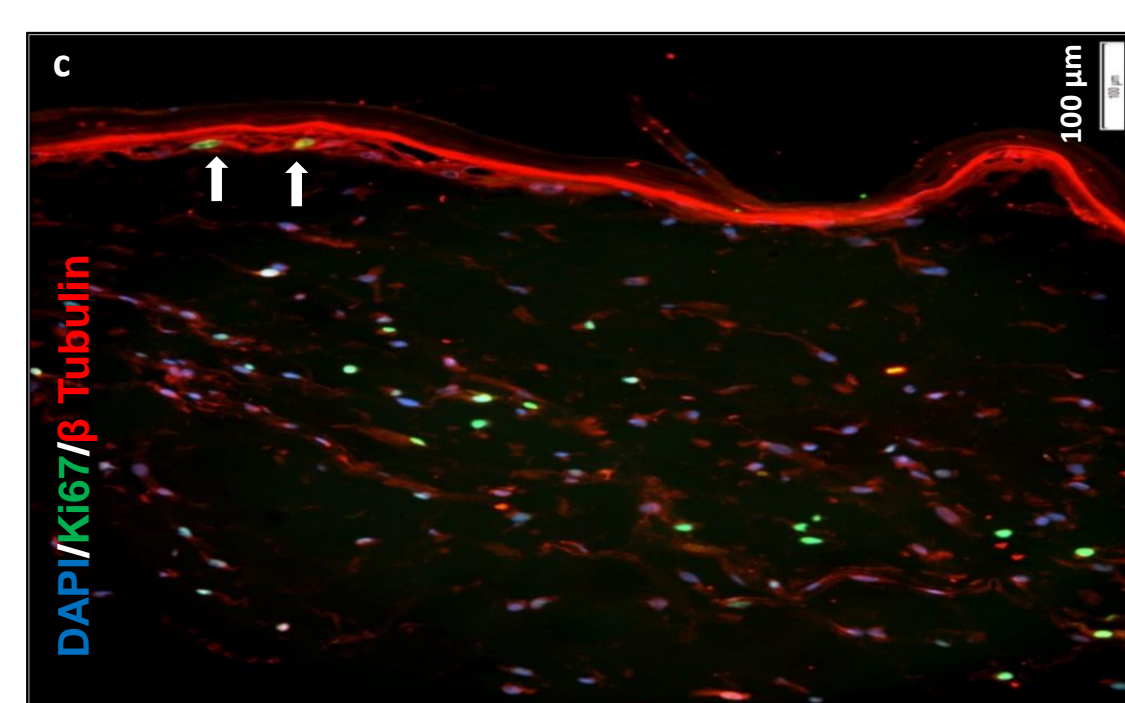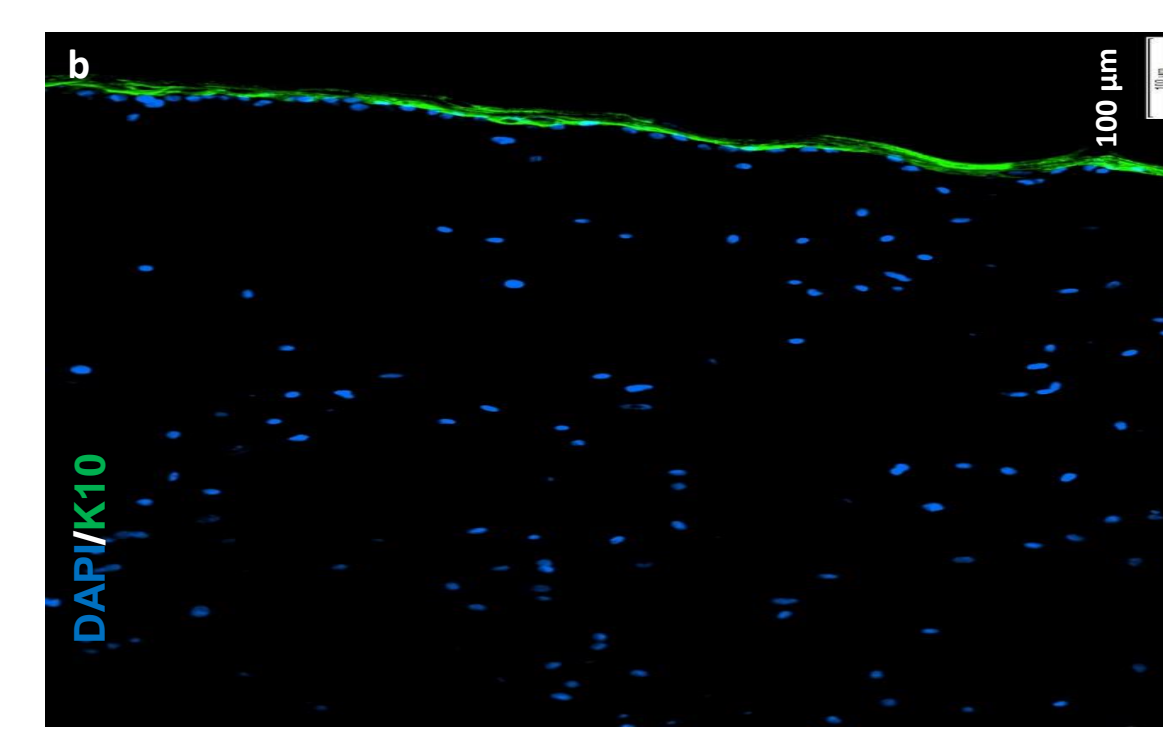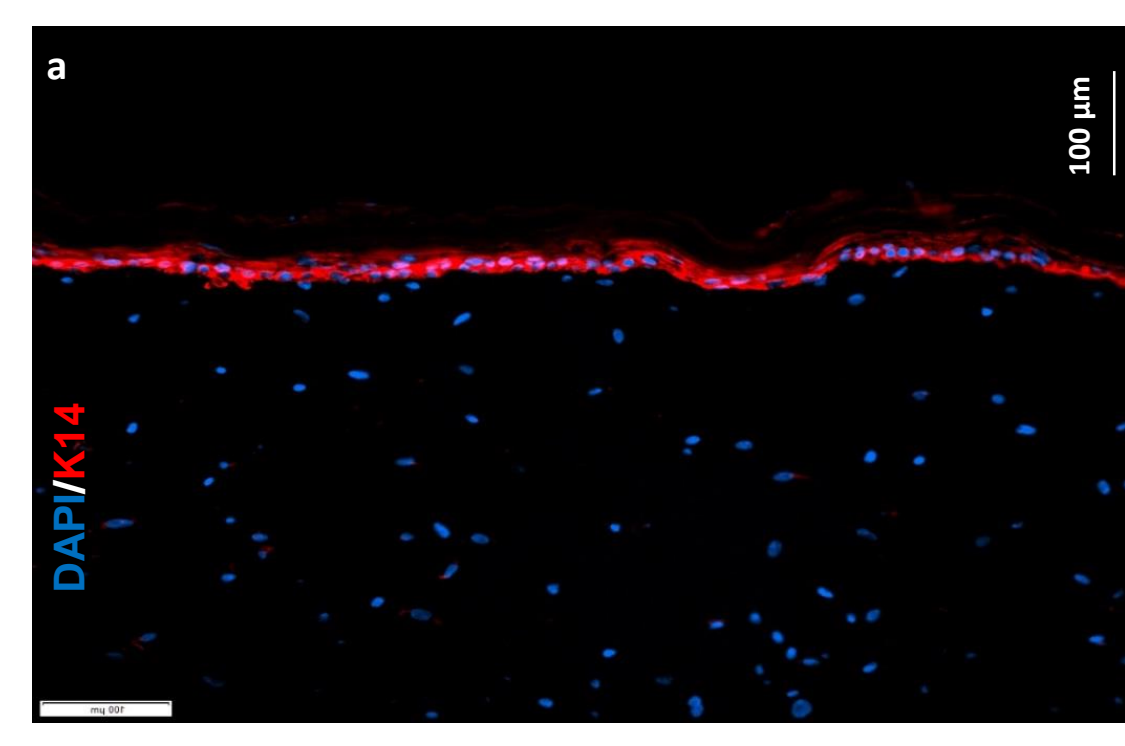


IHC verification of the Human Skin Equivalent (HSE):

- (a) Keratin 14 (K14) red;
- (b) Keratin 10 (K10) green;
- (c) Ki67 (arrows) green;
- (d) Laminin 5 (Lam5) green;

Cell nuclei are shown in blue by using DAPI staining. Laminin 5 is used as a marker of dermal-epidermal junction (DEJ) and appeared as a thin line.

Histological Analysis of the HSE

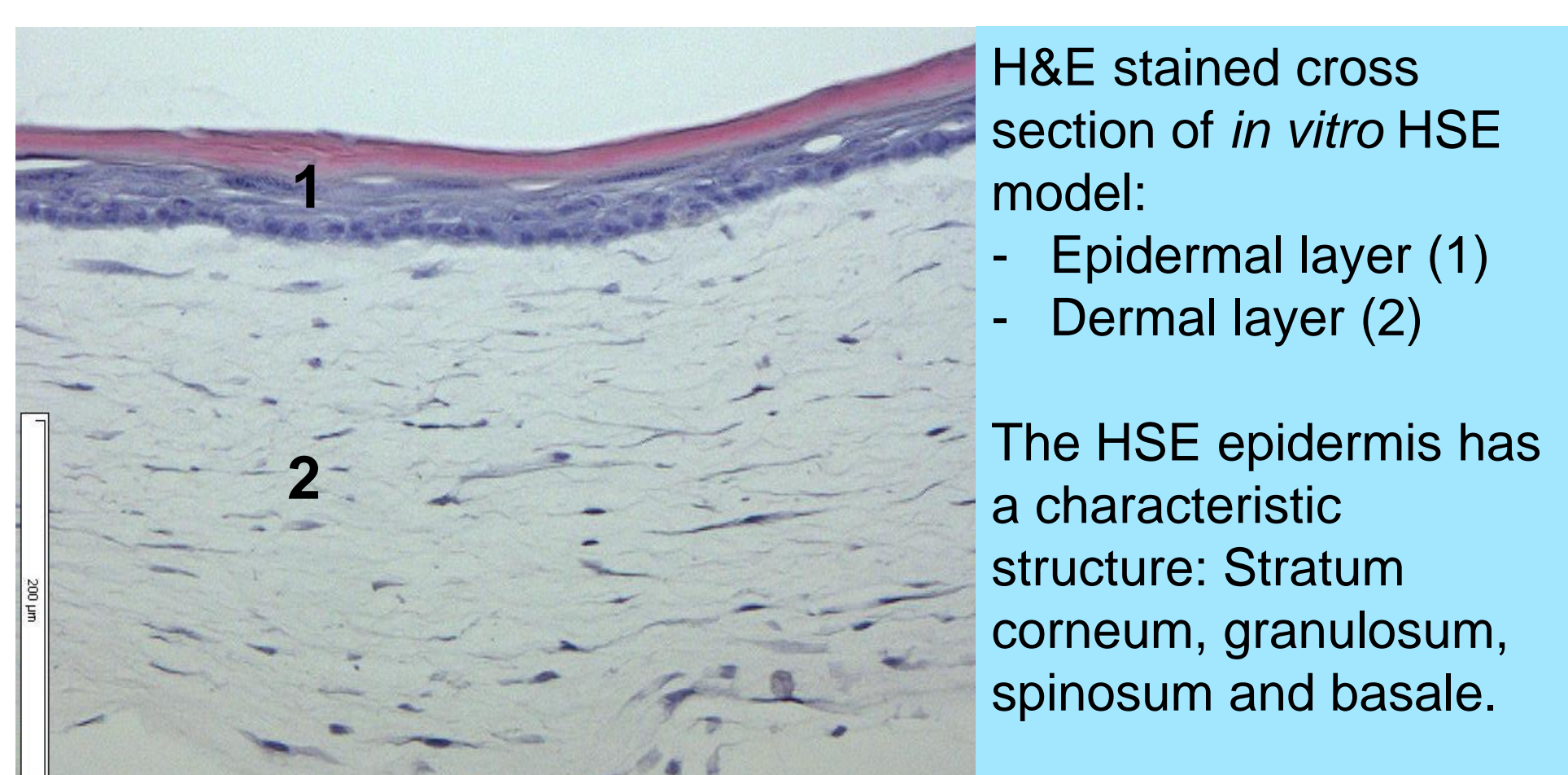


H&E stained cross section of *in vitro* HSE model:

- Epidermal layer (1)
- Dermal layer (2)

The HSE epidermis has a characteristic structure: Stratum corneum, granulosum, spinosum and basale.


Ultrastructure Analysis

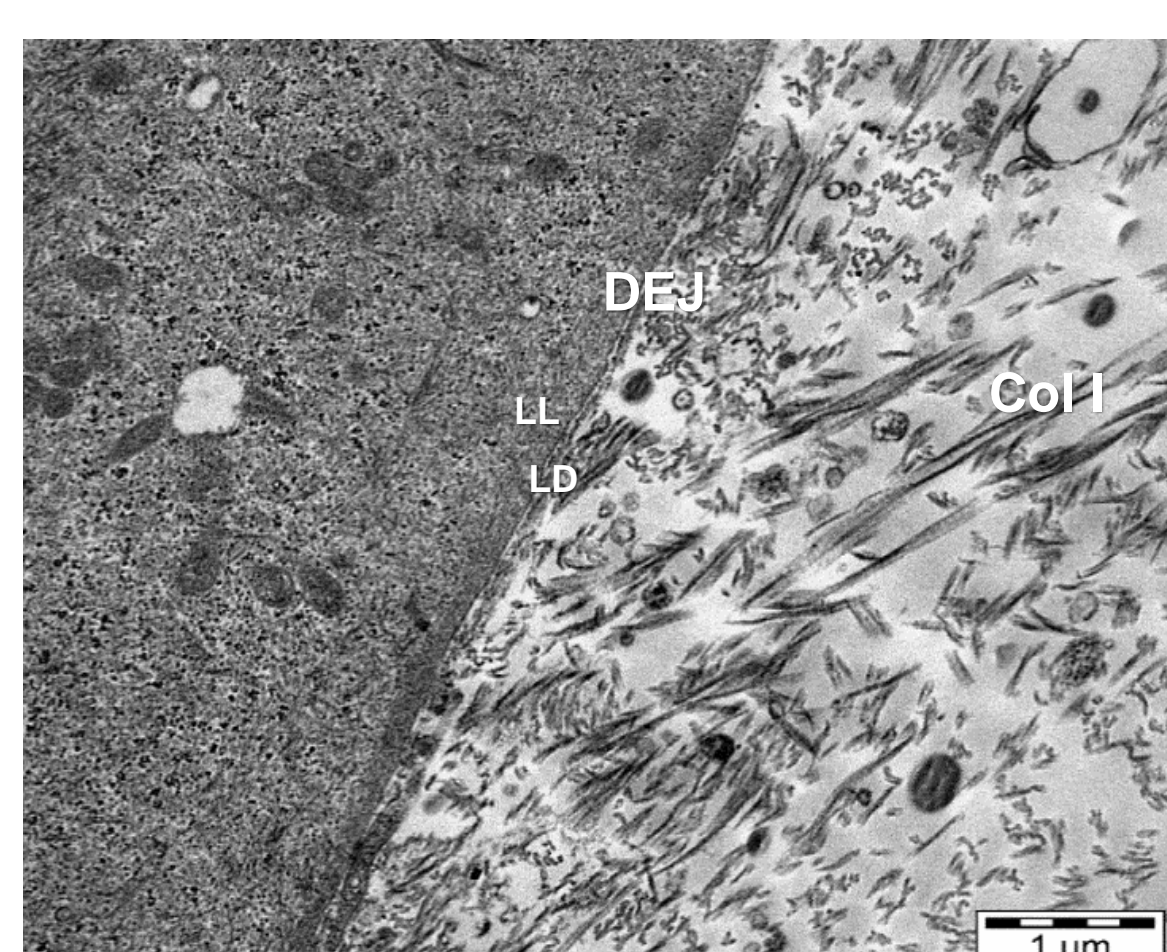

TEM image:

- Collagen-I fibres (Col I)
- Epidermal-Dermal Junction (DEJ).
- DEJ presents lamina lucida (LL) and lamina densa (LD).

S. aureus-colonized HSE

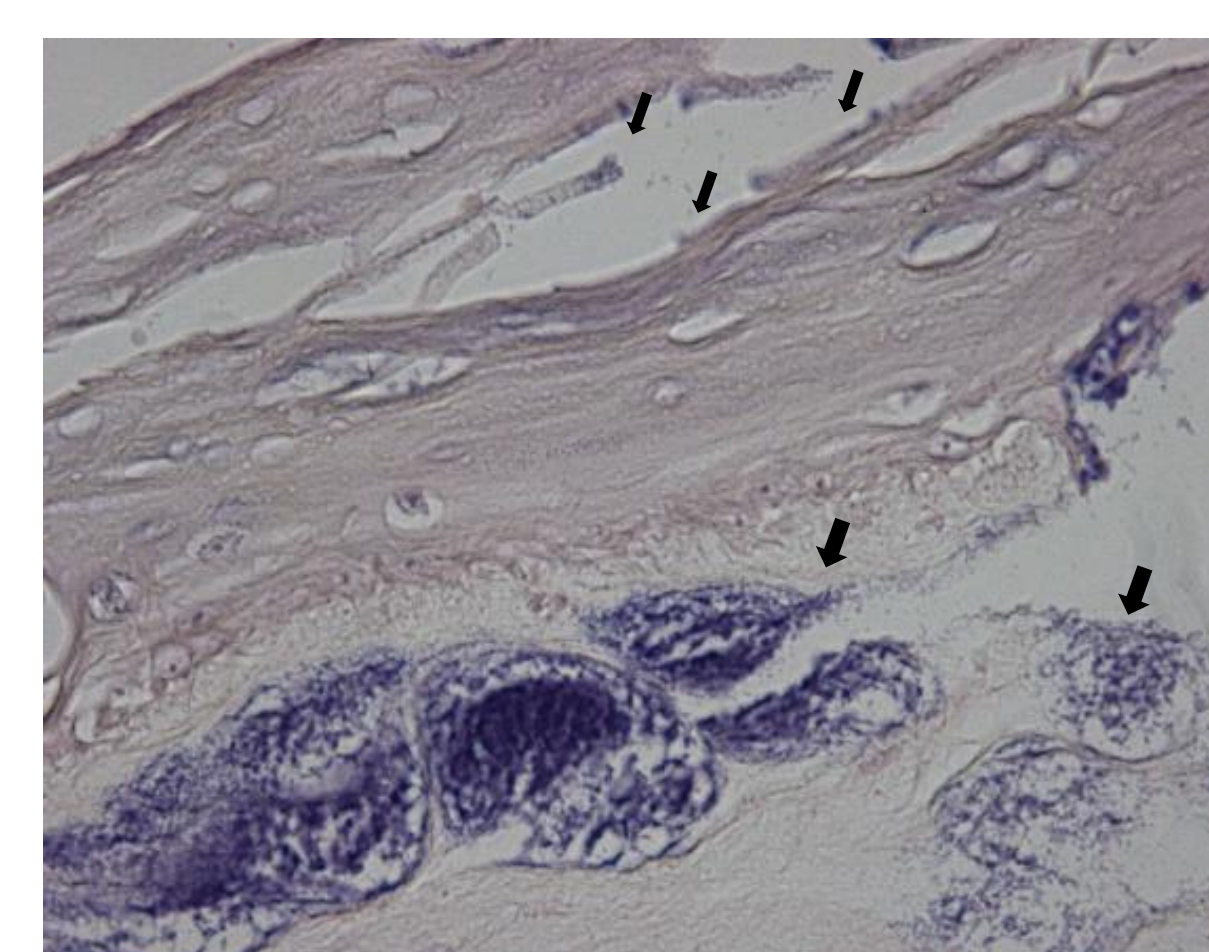


Inoculated bacteria adhere to the dermal surface, colonize, and replicate to make large structures of biofilm.

Big arrows: Bacteria located within a biofilm matrix inside the dermis.


Small arrows: Bacteria surrounding keratinocytes in epidermis.

Cytocompatibility Analysis

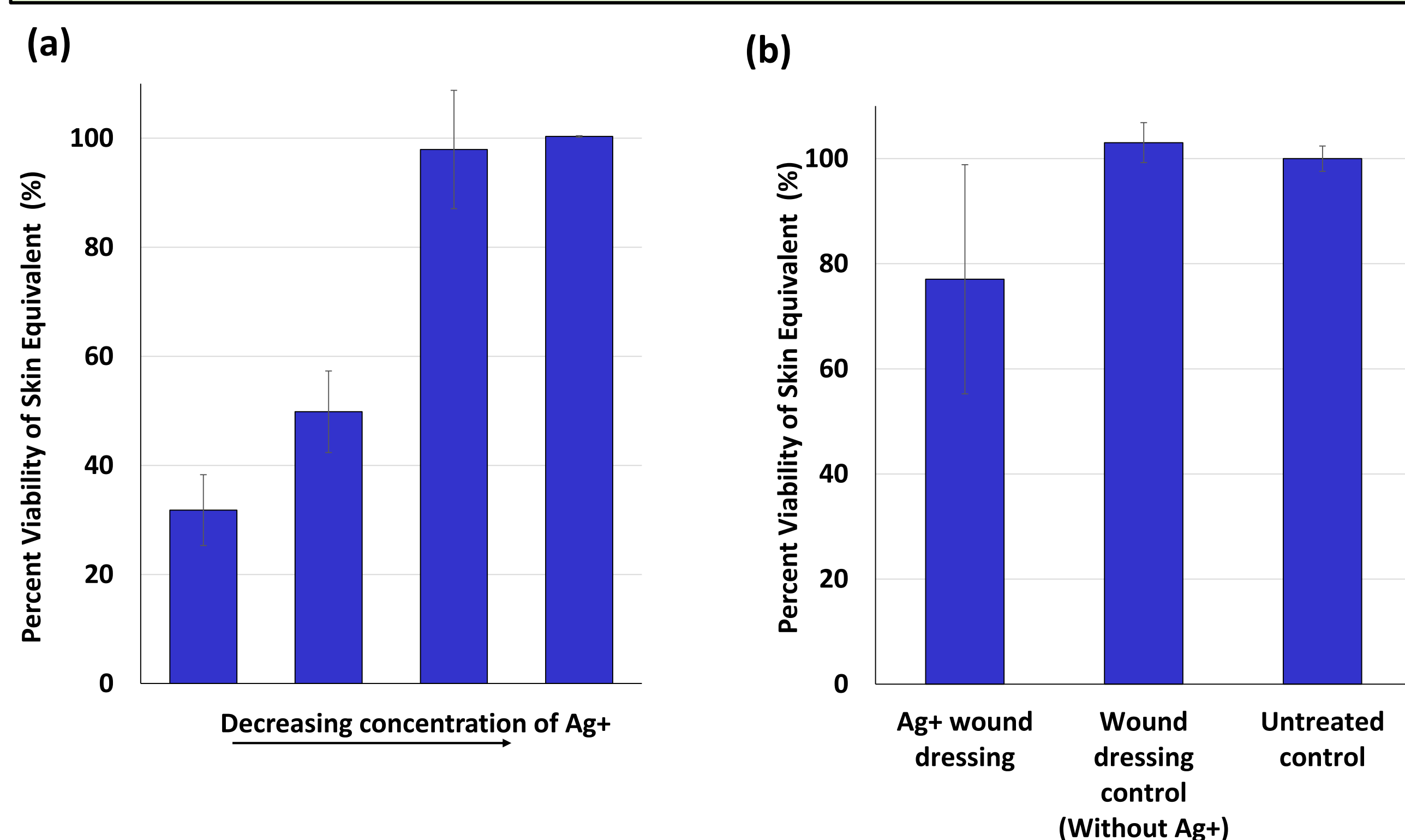

Cell viability measuring in the 3D system. The 3D skin model was exposed to a range of silver ion concentrations (Ag⁺) for a period of 24 hours. A commercially available Ag⁺ releasing wound dressing served the purpose of a model material and was tested in a 3D system along with its control material (without Ag⁺).

Antibacterial Analysis

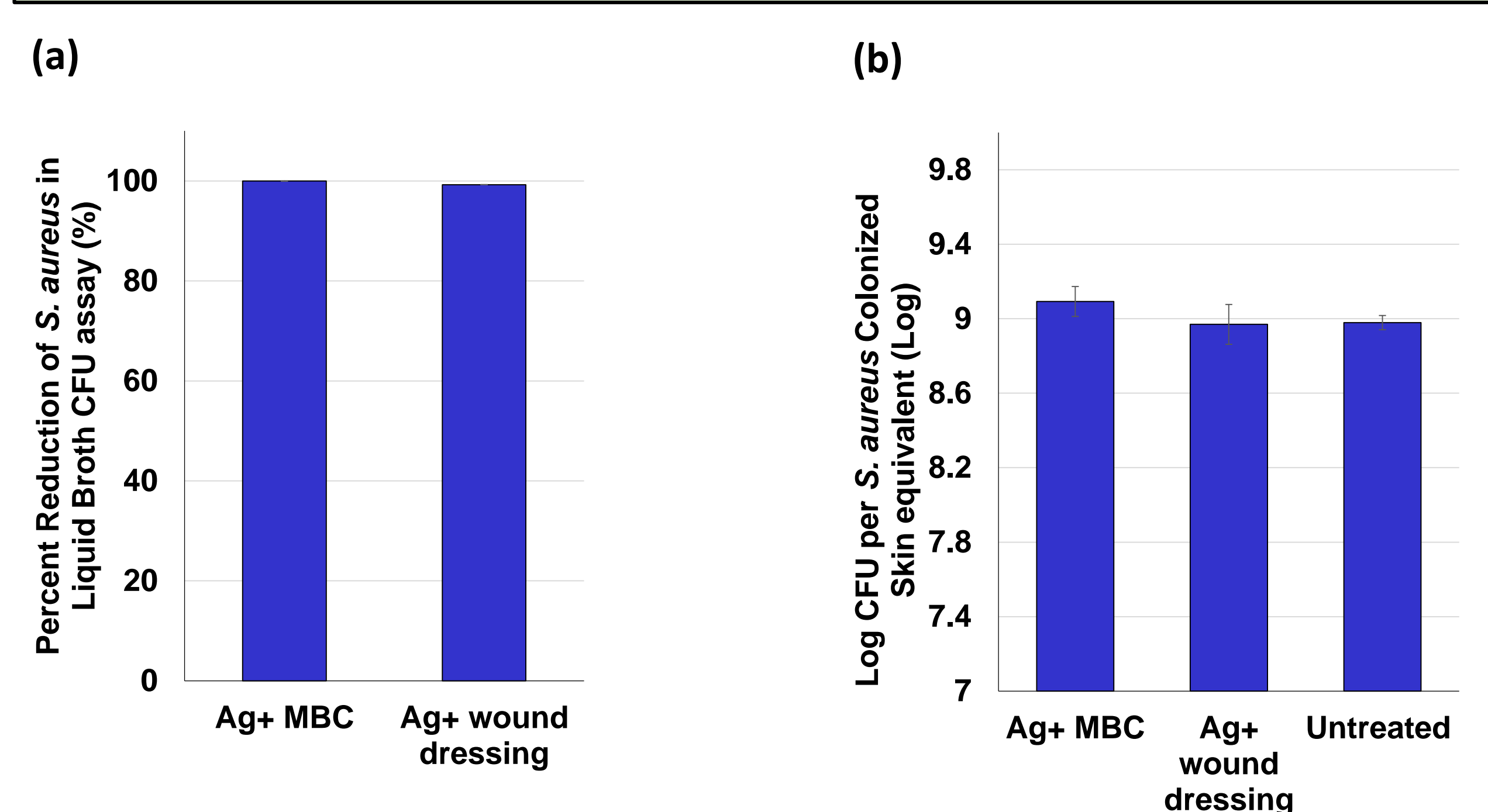

The graph demonstrates the treatment of infected skin equivalents with a commercially available Ag⁺ releasing wound dressing. Skin equivalents were infected with *S. aureus* and thereafter, Ag⁺ releasing wound dressing or Ag⁺ in PBS was applied onto the skin equivalents.

Conclusion: Development of colonized human skin equivalent (c-HSE); Risk assessment platform for cytocompatibility evaluation; Efficacy assessment of antibacterial materials; Comparison of 2D vs. 3D systems; Understanding "Host-Pathogen Interaction"; Development of complex skin models.